

ADVANCED

# MINISPIR<sup>TM</sup>

Handheld, PC-based Spirometer

Real-Time Flow/Volume and Volume/Time curves on your PC for a comprehensive Spirometry.



# MAIN features



## REAL-TIME TEST

Spirometry: FVC, VC, IVC, MVV, PRE/POST  
Bronchodilator comparison



## PLUG AND PLAY

Power via USB, no internal memory, no display, no maintenance, carrying case included



## COMPLIANCE ATS/ERS 2019

And other Standards including ISO 26782 (for Spirometry), ISO 23747 (for PEF), and more.  
CE0476, FDA 510 (k)



## SPIROMETRY PARAMETERS

Spirometry: FVC, FEV1, FEV1/FVC%, FEV3, FEV3/FVC%, FEV6, FEV1/FEV6%, PEF, FEF25, FEF50, FEF75, FEF2575, FET, ELA, EVOL, FIVC, FIV1, PIF, FIV1/FIVC%, PIF, IRV, VC, IVC, IC, ERV, FEV1/VC%, VT, VE, Rf, tI, tE, ti/tTOT, VT/tI, MVV



## PC CONNECTION VIA USB

Real-time test on PC screen, connect with your EHR/EMR, print Medical Report and more



# DISTINCTIVE features



## PREDICTED SETS & VALUES

Large Selection, including comparison %Pred, Z-score and LLN. Include GLI equations



## GENERAL PRACTICE

Easy-to-Use, real time spirometry curve and complete test results available in PC-mode



## EHR/EMR CONNECTIVITY

Via PC, integration with patient database on your EHR/EMR (in HL7, GDT)



## COVID-19 PREVENTION

Complete Disposable Set with Antiviral filter available, to reduce risk of cross-contamination



# Always INCLUDED

- Carrying case
- Noseclip
- PC Software license



# Compatible SOFTWARE

## winspiroPRO

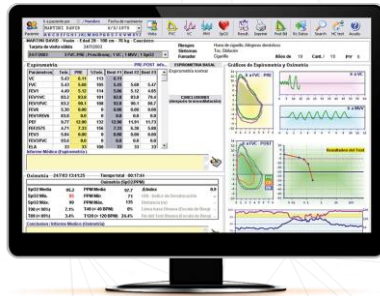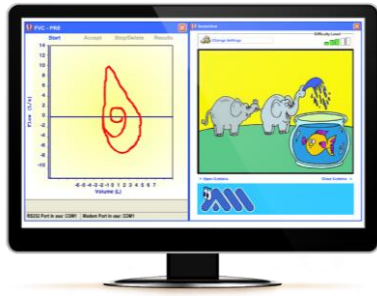

Pediatric Incentive (PATENTED) to improve patient compliance during the test.

Acceptability Messages, Test interpretation and Quality Control Grade according to the latest Spirometry Standards

### MAIN FEATURES

Windows-based solution for Spirometry, Oximetry and Telemedicine.


Wide range of predicted sets and values, including GLI Predicted sets, LLN and Z-score.

Embedded EHR/EMR connectivity.

NET VERSION available, share one database between different PC workstations.

### MEDICAL REPORT

Specialized and customizable printout

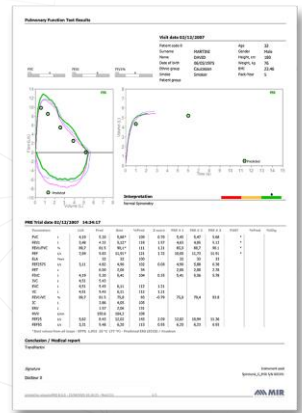

## spiro Connect



### MAIN FEATURES

Windows-based solution, direct integration with your EHR/EMR.

Real time test include Spirometry and Oximetry Standardized communication in HL7 or Exchange Protocol.

Select patient info directly from your own EHR/EMR

Spirometry test: FVC-Pre, FVC-Post, VC-Pre  
Oximetry Test: SpO2 (%), Pulse (BPM)

## GO-TO-MARKET TOOLKIT

Software Development Kit available for System Integrators and App Developers.  
OEM service available for Spirometry and Oximetry.



Learn more about available SDK and OEM



# Compatible TURBINES

flowMIR™ Disposable Turbine



Reusable Turbine



Mouthpiece

Included Disposable

Required, Not Included

Turbine Disinfection

Not required

Required

Turbine Calibration

Not required

Required

Packaging

Individually sealed: 60 or 10 units / box

1 unit in Carton box

Antiviral Filter

Available Disposable

Required Disposable

PLAY VIDEO



SCIENTIFIC PUBLICATIONS

